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Device for improvement
of own defenses and
correction of pure
electromagnetic radiation
emitted by natural biosystems

DeVita
Energy

Intended Use

«DeVita Energy» is a generator of low-level electromagnetic radiation identical to emissions of natural organic systems.

The device is designed to improve the body's own defenses and correct pure electromagnetic radiation emitted by natural biosystems. The principle of corrective action is based on the usage of R.Rife, R.Voll, and F.Morrel frequencies.

Caution! This device can be used by children older than 8 years old and by persons with limited capabilities such as physical, sensory or other mental health issues. In addition, it can be used by persons with a lack of experience and knowledge provided that the user will be under supervision or the user is provided instructions concerning use of the device and understands the associated risk. Children are forbidden to play with the device. Cleaning and servicing the device should not be carried out by children without supervision.

Content of Delivery

«DeVita Energy».....	1 pce
Certificate/Instruction for Use.....	1 pce
Package.....	1 pce

Specification

1. Complexes, pcs	1
2. Frequency range, Hz.....	0,1÷15000
3. Frequency setting error, max, %.....	±0.003
4. Run time, min, hours.....	20
5. Keyboard type.....	touch-sensitive keyboard
6. Display type.....	TFT 128x160
7. Online update for:	
- device control program	available
- complex programs	available

8. Power source:	
- Li-Pol 3.7V Battery, device-integrated, pce	1
9. Battery maximum charge time, hours	6
10. Dimensions, max, mm.....	102x44x15
11. Device weight, max, kg.....	0.1
12. Mean life time, min, years.....	5

Safety

The device is:

- safe to use
- non toxic
- non flammable
- non explosive
- not a source of radiation hazardous to health.

When charging, III class devices use SELV (Safety Extra Low Voltage) network adapters only.

Environmental Protection

The integrated battery contains environmental contaminants. After expiration of the devices useful life, the battery must be removed before disposal. The battery must be disposed at special stations. 

Battery Removal

1. Remove three screws on the back panel of the device body.
2. Disconnect sealed battery terminals from the board.
3. Remove the battery with care.

Instructions for Service, Storage and Transportation

The unpacked device must be kept indoors at temperatures ranging from 10 - 35°C and a maximum relative humidity of 80%.

The device must be kept away from direct moisture and exposure to chemical agents

Do not disassemble the device.


Prevent the device from falling.

Prevent objects from falling onto the device.

We recommend keeping the package during the entire warranty period.

Device Description

1. ON/OFF Button
2. Display
3. Charging socket
4. Direction key
5. Option selection key
6. Return key
7. Program execution LED

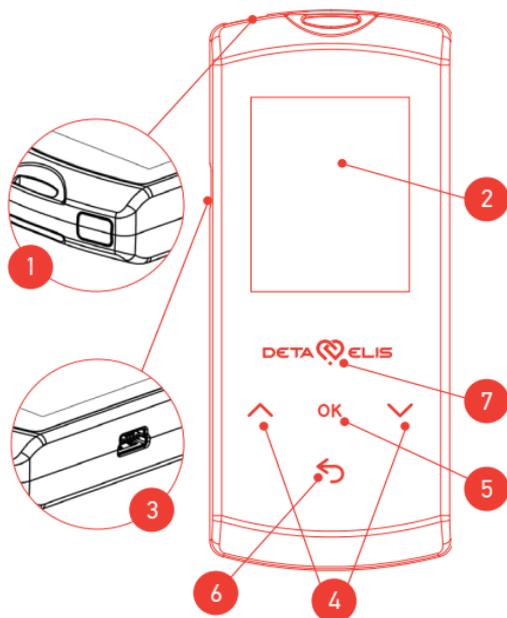Operation

Press button¹ to switch on the device. The display will remain lit for 30 seconds from the moment of pressing the last key. The display can be turned off earlier by pressing button¹ for the second time. If the display turns off before you enter the required commands, press button¹ again to turn on the display.

After pressing button¹, the start up screen¹ appears on the display and in 2 secs the device will automatically enter the main menu² mode.

Working with Main Menu

"Setting" Submenu

Select "Setting" item with the direction keys . Press  key to select options.

"Setting" submenu³ will appear on the display.

• Language Setting

Attention! The Russian language setting is set by default.

If there is no need to change the language, skip this item.

Select "Language" item with the direction keys . Press  key to select options. "Language" submenu⁴ will appear on the display.

Select the required language with the direction keys . Press  key to select options.

To return to "Settings" submenu, press  key.

• Audio Tuning

Select "Audio" item with the direction keys . Press  key to select options. "Audio" sub-menu⁵ will appear on the display.

Select "On" or "Off" item with the direction keys . Press  key to select options. Selecting "Off" item will switch off the key audio feedback.

To return to "Settings" submenu, press  key.

• Time Setting

Select "Time" item with the direction keys . Press  key to select options. "Time Setting" sub-menu⁶ will appear on the display.

Select "Hours", "Minutes" or "Format" window with the direction keys .

Set the required hours, minutes or format with direction keys .

To return to "Settings" submenu, press  key.

• About Device


Select "About device" menu item with the direction keys . Press  key to select options. The number of your device and program version will appear on the display.

• Program Selection


Select "Program" item² with direction keys . Press  key to select options. The blinking of the red light on the device panel under the logo indicates the program will launch. The "current program execution" submenu⁸ will appear on the display. There is a timer located on the upper part of the display indicating the remaining program execution time.

1

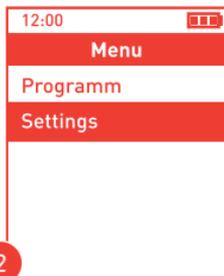2

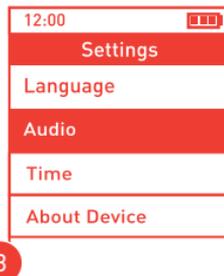3

4

• Information about Program Execution

If any information about program execution is required, press button¹. The program execution sub-menu will appear on the display.

• Program Interruption

If the program must be interrupted, select the "Yes" item of from the "Program execution" sub-menu by using the direction keys   and press  key to select options.

• Use of Device in Program Working Mode

Immediately after program launch, place the device on the body. For example, put the device into your pocket.

Device Charge

• Charge Indicator

The Charge Indicator is white-lighted, two  or three  internal indicator sectors are active – the battery is charged.

The Charge Indicator is red-lighted, only one internal indicator sectors  is active – the battery is low.

The Charge Indicator is red-lighted, no internal indicator sector is active, "Charge the battery" is on the display – the battery is in need of charge⁹.

• Battery Charging

If the battery charge is low or the battery is in need of charge, connect the provided mini-USB connection cable to the USB port of the computer. The charge indicator will indicate charging is taking place.

Manufacturer's Warranty

The manufacturer guarantees compliance of the DeVita Energy device to the technical specifications

provided to user for operation, transportation and storage.

The covered warranty period of the device is 24 months starting from the date of sale through the retail chain. If the device certificate does not have a sale date or a seller stamp, the warranty period starts from the date of the production of the device.

The user has the right for the free repair within the warranty period. All repairs under warranty shall be performed at the manufacturer's facility. The defective device shall be brought to the repair station at the buyer's expense.

In case of device body impairment and/or absence of warranty certificate, the claims are not valid and the warranty repair will not be performed.





The warranty does not cover the following defects:

- defects due to improper operation
- defects due to natural disasters
- in case of broken seals
- in case of exterior defects (cracks, cleavages, and others).

The buyer has the right to have the defective device substituted for the new one in the following cases:

- the device was under repair three times during the warranty period;
- the device is beyond repair.

Attention! The manufacture reserves the right, without any prior notice, to make modifications to the device design and software without compromising current parameters.

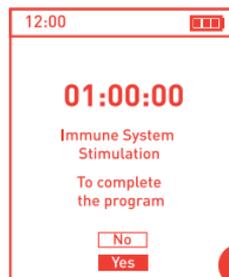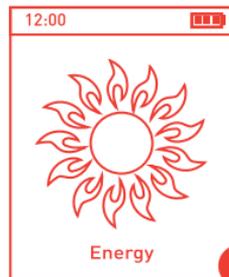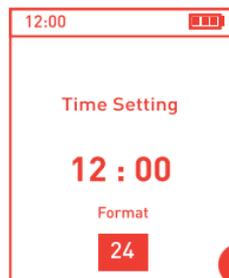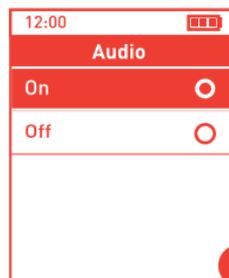Certificate of Acceptance

"DeVita Energy"

serial No. _____,

is manufactured and accepted in accordance with the obligatory state standards, effective technical documentation, and classified as fit for operation.

Technical Control Mark:

Seal here _____
(personal signature) (print full name)

« _____ » _____ 201__ г.
(issue date)

For the first time, the notion of the existence in the human body of energy and health's dependence on body circulation was stated by ancient medical healers in China, India, and Egypt. Today, it is scientifically proven that each human organ radiates certain frequencies. Because of the influence of adverse internal factors, modern life stressors, and a polluted environment, frequencies often deviate from their natural parameters. As a result, depression, apathy, and chronic sleeplessness occur.

"DeVita Energy" devices radiate harmonious electromagnetic frequencies which balance biopower indicators for each person, stimulate functions for body regulation, based on the resonant phenomena and protect from biological and geopathogenic factors. "DeVita Energy" is soft and safely influences a human body, without breaking its energy balance as a whole or on a local level.

This way of understanding the body is based on the work of Franz Morell, the author of bioresonant therapy methods. A correctly created program that regulates electromagnetic frequencies restores and strengthens physiological indicators, removes psycho-emotional stress from a person, and increases a person's energy level providing a person with a good mood and fine health.

Attention! As a safety precaution, this device should not be used by pregnant women in the first trimester, people who have received organ transplants, and also for people with serious heart conditions.

Indications for use of "DeVita Energy" device

It is recommended that people use the device who suffer from the following conditions:

- Physical and emotional fatigue
- Chronic fatigue syndrome
- Large stress load
- Depression, Stress
- Dizziness, weakness, drowsiness

Effects after using the "DeVita Energy" device:

- Restores and strengthens physiological indicators
- Reduces psycho emotional stress

- Increases energy level
- Stimulates functions for body regulation
- Regulates drainage function of the body
- Strengthens biological and geopathogenic factors
- Rejuvenates

Technique: Use the device 1-2 times a day in the morning or day time. Do not use the device at night or more than 2 times a day.

Conclusion

Wellness concept – healthy food, moderate physical stress burden, fresh air, abstaining from unhealthy habits – in many cases, this reduces risk of most diseases. However, in the modern world, it is not always possible to adhere to this concept completely. In particular, it is common for large cities to be filled with sources of stress, a fast-paced way of life, and an unfavorable environmental situation.

"DeVita Energy" devices are needed to help the body restore morning energy, increase immunity and to support overall life vitality on a high level.

The more energy a person has, the better his state of health. This will provide for increased potential and opportunities for a person to create the specific lifestyle that he/she wants.

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DETA-ELIS-EUROPA GmbH Justus-Liebig-Str. 2
DETA-ELIS-Europa GmbH